Shower seat

LinidoSolutions™

Manual - English
Manual - Svenska
Manual - Norsk
Manual - Dansk
Käyttöohje - Suomi
Gebrauchsanweisung - Deutsch
Handleiding - Nederlands
Manual - Français
Manual - Italiano
Manual – Español



LI2202.2006-xx




LI2201.2006-xx



LI2203.2006-xx

Art.no.:
LI2201.2006-xx
LI2202.2006-xx
LI2203.2006-xx

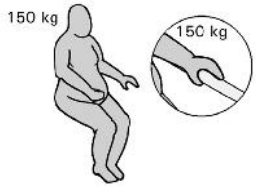

Max: 150 kg/330 lbs

handicare

Manual no: LI09509.461 140224



EN Visual inspection
Inspect the product regularly. Check to ensure that material is free from damage.

SE Visuell inspektion
Gör regelbundna funktionskontroller av produkten. Kontrollera att materialet är helt fritt från skador.

NO Visuell inspeksjon
Gjør regelmessige funksjonskontroller av produktet. Kontroller at materialet er helt og uten skader.

DK Visuel inspektion
Foretag regelmæssige funktionskontroller af sengebordet. Kontroller, at materiale er helt og ikke udviser tegn på beskadigelse.

FI Silmämääräinen tarkastus
Tarkasta tuotteen toimivuus säännöllisesti. Tarkasta, ettei materiaalissa ole mitään vaurioita.

DE Sichtprüfung
Unterziehen Sie das Hilfsmittel regelmäßigen Funktionskontrollen. Vergewissern Sie sich, daß das Material unbeschädigt ist.

NL Visuele controle
Inspecteer het product regelmatig. Controleer op mogelijke beschadigingen van het materiaal.

FR Inspection visuelle
Effectuer des inspections régulières du fonctionnement du produit. Vérifier que le matériau est entièrement exempt de dommages.

IT Ispezione visiva
Testare il funzionamento del prodotto regolarmente. Verificare che il materiale non sia danneggiato.

ES Inspección visual
Realice periódicamente un control de funcionamiento del producto. Asegúrese de que el material no presente ningún tipo de desperfecto.

EN Before use
Ensure that the product is correctly mounted/installed.

SE Före användning
Kontrollera att produkten är rätt monterad/installerad.

NO Før bruk
Kontroller at produktet er korrekt installert.

DK Før brug
Kontroller, at produktet er korrekt installeret.

FI Ennen käyttöä
Tarkasta, että tuote on asennettu oikein.

DE Vor der Benutzung
Prüfen Sie das Produkt auf korrekte Montage.

NL Voor gebruik
Zorg ervoor dat het product juist is gemonteerd.

FR Avant l'utilisation
Vérifier que le produit est correctement installé.

IT Prima dell'uso
Verificare che il prodotto sia installato correttamente.

ES Antes de su uso
Compruebe que el producto esté correctamente instalado.



EN Always read the manual
Keep the manual where it is easily accessible for users of the product.

SE Läs alltid manualen
Förvara manualen tillgänglig för användare av produkten.

NO Les alltid manualen
Oppbevar manualen tilgjengelig for brukere av produktet.

DK Læs altid manualen
Opbevar manualen tilgængelig for brugere af produktet.

FI Lue aina käyttöohje
Säilytä käyttöohje tuotteen käyttäjien saatavilla.

DE Lesen Sie stets die gebrauchsanweisung
Verwahren Sie die Anleitung so, dass sie der Anwender des Produktes rasch zur Hand hat.


NL Lees altijd de handleiding
Bewaar de handleiding zo dat deze voor gebruikers van het product altijd beschikbaar is.

FR Toujours lire le manuel
Conserver le manuel pour qu'il soit accessible aux utilisateurs du produit.


IT Leggere sempre il manuale
Tenere il manuale delle istruzioni in un luogo facilmente accessibile agli operatori che utilizzano il prodotto.


ES Lea siempre el manual
Guarde el manual en un lugar accesible para el usuario del producto.


Using/Användning/Bruk/Brug/Käyttö/Anwendung/Gebruik/Utilisation/Utilizzo/Empleo

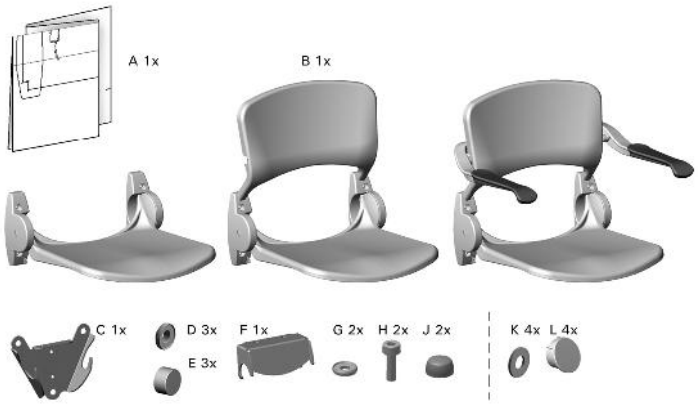



GB: Use proper fixing materials for the wall and resistant to corrosion.
SE: Använd skruvar av rostfritt material. Skruvarna ska vara anpassade till materialet i väggen.
NO: Bruk rustfrie skuer. Skruene skal være tilpasset materialet i veggen.
DK: Använd materialer til væggen, som er korrosions frie.
FI: Utiliser des chevilles adaptées au mur et des fixations inoxydables.
DE: Benutzen Sie Schrauben aus rostfreiem Material und die für die Wand geeignet sind.
NL: Gebruik bevestigingsmaterialen geschikt voor de muur en bestand tegen corrosie.
FR: Utiliser des chevilles adaptées au mur et des fixations inoxydables.
IT: Usare materiali di fissaggio a parete adatti al tipo di muro e resistenti alla corrosione.
ES: Usar tornillos apropiados para cada tipo de pared y que sean anticorrosivos (inox).

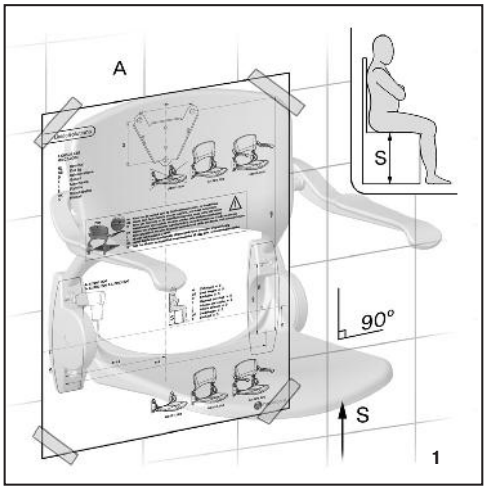


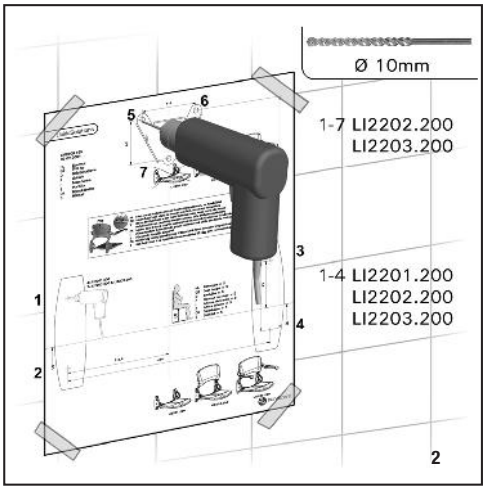


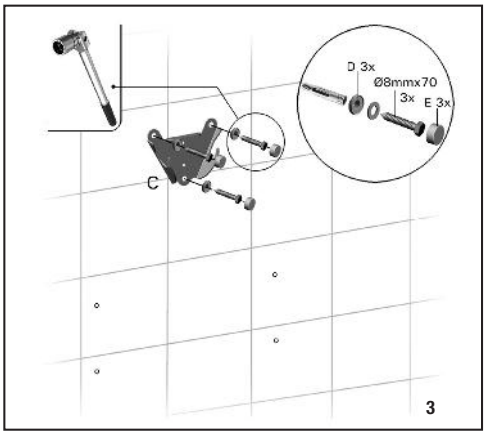

D 3x
Ø8mmx70
3x
E 3x
C
3

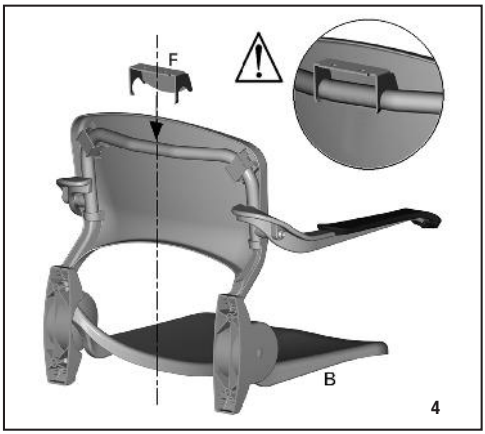
F
B
4

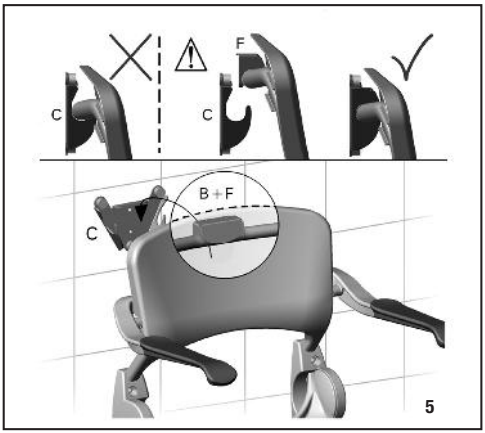

C
F
C
B+F
C
5

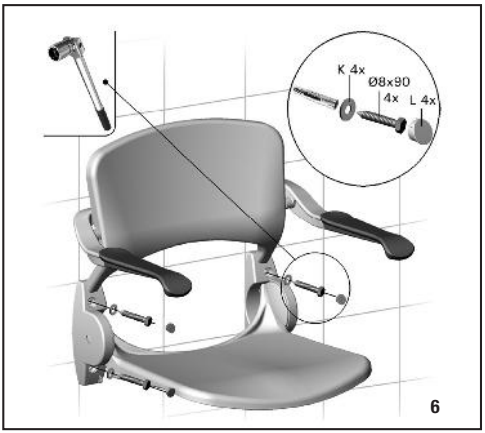

K 4x
Ø8x90
4x
L 4x
6

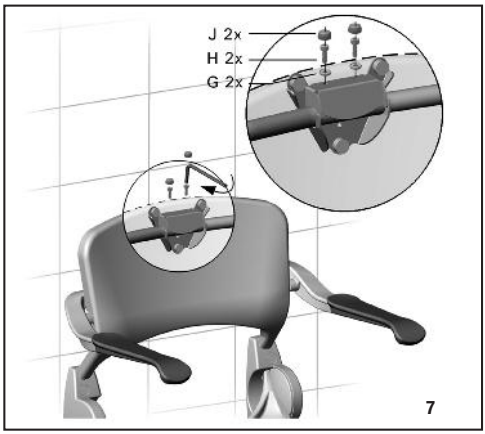
J 2x
H 2x
G 2x
7

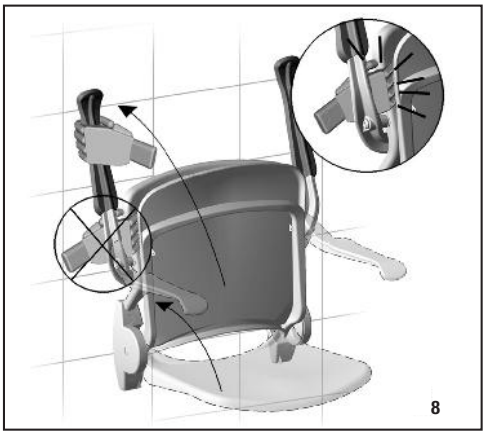
8



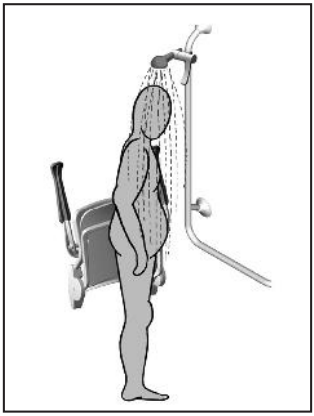

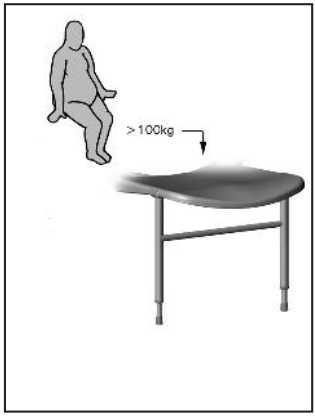

Maintenance/Rengöring/Vedlikeholdsråd/Vedligeholdelse/Hoito-ohje/Pflegeanleitung/ Onderhoud/Conseils d'entretien/Cura del prodotto/Ciudado del producto

EN Clean regularly with normal household detergent. Do not use aggressive cleaners. Handicare is not liable for damage caused by improper mounting/use.

SE För rengöring använd normalt och inte aggressivt rengöringsmedel. Handicare står ej ansvarig för skador som uppstår pga. felaktig montering/användning.

NO For rengjøring bruk nøytralt/normalt rengjøringsmiddel og ikke etsende midler. Handicare står ikke ansvarlig for skader som kan oppstå pga feilaktig motering/bruk.

DK Efter brug tjernes sæberester. Dar må ikke anvendes hårde rengøringsmidler. Handicare er ikke ansvarig for skader opstået pga. torkert montering/brug.

FI Puhdista säännöllisesti käyttäen kodin yleispesuaineita. Älä käytä voimakkaita liuottimia. Handicare ei ole vastuussa vahingosta, joka johtuu virheellisestä asennuksesta tai käytöstä.

DE Beseitige nach Benutzung Seifenreste. Kein agressieves Reinigungsmittel verwenden. Handicare ist nicht haftbar für Schäden die durch unsachgemässe Befestigung/Benutzung oder falsches Befestigungsmaterial entstanden sind.

NL Verwijder na gebruik zeepresten. Gebruik geen agressieve schoonmaakmiddelen. Handicare is niet aansprakelijk voor schade als gevolg van onjuiste bevestiging/gebruik.

FR Laver régulièrement avec un détergent ménager. Ne pas utiliser de détergent agressif. Handicare n'est pas responsable des dommages causés par un montage/utilization incorrect.

IT Pulire periodicamente, utilizzando detergenti di uso corrente per la casa. Non usar prodotti abrasivi e corrosivo. Handicare non è responsabile per i danni causati dal montaggio/uso impropio.


ES Limpiar regularmente con cualquir detergente no corrosivo. No usar leija ya que es muy corrosivo. Handicare no se hace responsable de cualquir deterioro causado por un montaje/uso incorrecto.

Accessories/Tillbehör/Tilbehør/Lisävarusteet/Zubehör/Accessoires/Accessori/Accesorios

Art. no:	Product namne:	Description:
LI2806.0000-00	Wall mounting set	For wall mounting of all variants of Ergonomic shower seat without backrest.
LI2805.0000-00	Wall mounting set	For wall mounting of all variants of Ergonomic shower seat with backrest.
LI0122.5084-02	Backing plates set	For secure fixing to a weaker wall, for all variants of Ergonomic shower seat.
LI09506.400	Cushioned cover	Soft seat cover, for all variants of Ergonomic shower seat
LI2216.0014-xx	Click-on frame	For wall mounting of all variants of Ergonomic shower seat with backrest, for hanging on Ergogrip 600/700/800/1600 and all variants of Ergogrip corner.
LI2216.0024-xx	Click-on frame	Height adjustable, for all variants of Ergonomic shower seat with backrest, for hanging on Ergogrip 600/700/800/1600 and all variants of Ergogrip corner.
LI2214.201x-xx	Floor support for Ergonomic shower seats	For all variants of Ergonomic shower seats, for users weighing more than 100kg.
LI2215.1014-xx	Height adjustment system	For wall mounting of all variants of Ergonomic shower seat without backrest.
LI2215.1024-xx	Height adjustment system	For wall mounting of all variants of Ergonomic shower seat with backrest.
LI2603.3xxx-xx	Hinged arm support	For side support in combination with all variants of Ergonomic shower seat without armrests, for larger users weighing more than 125 kg.
LI2801.0000.02	Cover set	Spare part, for Ergonomic shower seat LI2201.2006-xx, white.
LI2803.0000-02	Cover set	Spare part, for Ergonomic shower seat LI2202.2006-xx and LI2203.2006-xx, white.


Dimensions/Måttbeskrivning/Målbeskrivelse/Mitat/Maßbeschreibung/Afmetingen/ Dimensioni/Dimensiones

Ergonomic shower seat | Art. no. LI2201.2006-xx

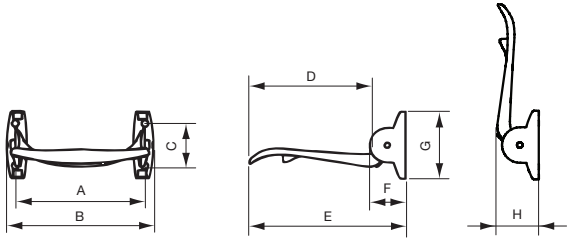

A	425 mm	16,7 inch
B	485 mm	19 inch
C	140 mm	5,5 inch
D	405 mm	15,9 inch
E	515 mm	20,2 inch
F	120 mm	4,7 inch
G	220 mm	8,6 inch
H	140 mm	5,5 inch


Ergonomic shower seat, with backrest | Art. no. LI2202.2006-xx

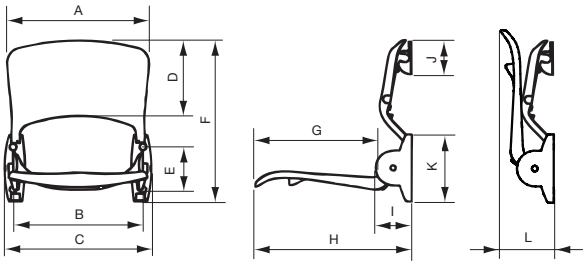

A	465 mm	18,3 inch
B	425 mm	16,7 inch
C	485 mm	19 inch
D	245 mm	9,6 inch
E	140 mm	5,5 inch
F	530 mm	20,8 inch
G	405 mm	15,9 inch
H	515 mm	20,2 inch
I	120 mm	4,7 inch
J	120 mm	4,7 inch
K	220 mm	8,6 inch
L	170 mm	6,7 inch

Ergonomic shower seat, with backrest and armrests

Art. no. LI2203.2006-xx

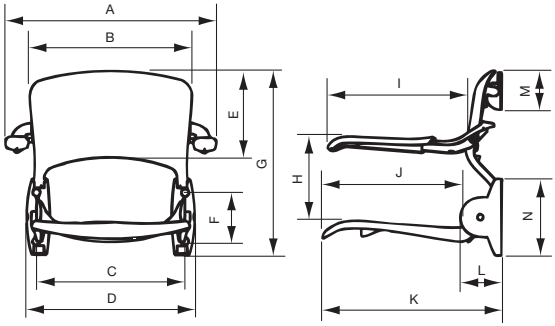

A	605 mm	23,8 inch
B	465 mm	18,3 inch
C	425 mm	16,7 inch
D	485 mm	19 inch
E	245 mm	9,6 inch
F	140 mm	5,5 inch
G	530 mm	20,8 inch
H	260 mm	10,2 inch
I	405 mm	15,9 inch
J	405 mm	15,9 inch
K	515 mm	20,2 inch
L	120 mm	4,7 inch
M	120 mm	4,7 inch
N	220 mm	8,6 inch

LinidoSolutions™

Handicare's LinidoSolutions™ product line is a wide range of high quality safety products designed to provide ergonomic support, comfort and increased independence for the user in various bathing, showering and toileting situations.

Contact your local distributor if you have any questions about the product and its use. See www.handicare.com for a complete list of distributors. Always make sure that you have the right version of the manual. The most recent editions of manuals are available for downloading from our website, www.handicare.com

MediCauce

MediCauce Soluciones Sanitarias SL
C/ Rio Jarama 132, 10-13
45007, Toledo. España
T- (+34) 925 337 467
www.medicauce.com
info@medicauce.com